CoCrMo-Sold 1

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This product contains cobalt.

Free from nickel, beryllium, lead, cadmium and iron.

Technical data

Solidus temperature	T _s	1110 °C / 2921 °F
Liquidus temperature	T _l	1162 °C / 3002 °F
Operating temperature (recommended)	T _a	1262 °C / 2309 °F

Dear customer,

Thank you for choosing a quality product from Dentaurum. It is essential to read these instructions carefully and adhere to them to ensure safe, efficient use and ensure that you and your patients gain full benefit. Instructions for use cannot describe every eventuality and possible application. In case of questions or ideas, please contact your local representative.

As our products are regularly upgraded, we recommend that you always carefully read the current Instructions for use supplied with the product and stored in the internet at www.dentaurum.com, even though you may frequently use the same product.

1. Manufacturer

Dentaurum GmbH & Co. KG | Turnstraße 31 | 75228 Ispringen | Germany

2. General device description

Dental solder for cobalt-based alloys.

3. Area of application

CoCrMo-Sold 1 is a dental solder for cobalt-based cast alloys. It is particularly suitable for soldering alloys (e.g. remanium® GM 380+, remanium® GM 800+, remanium® secura or remanium® 2000+) which are intended for the production of crowns and bridges and metal ceramics. It is also suitable for partial denture alloys based on cobalt.

4. Recommended flux

For soldering, rema® Flux 1 (REF 102-304-00) is recommended. It is a flux specially adapted to the dental solder.

5. Preparation

In preparation, a solder base made of an investment material for partial dentures, e.g. rema® Exakt F or rema® dynamic S, with optimally adapted compensation of the metal contraction (Fig. 1) should be used. Water is recommended as the mixing liquid. In general, the solder block should be made as small as possible and as stable as necessary. The soldering gap is ideally 0.2 mm and should not exceed 0.5 mm. The solder block should be dried a temperature of 200 °C / 392 °F for a period of 30 minutes.

Alternatively, instead of the soldering model, the bridge can be dotted „triangle-shaped“ with metal wires of the same type (Fig. 2).

6. Procedure

A solder burner with sufficient power (e.g. propane/oxygen) should be used. The solder joint is covered with the flux rema® Flux 1. The solder area should be covered by the reduced zone of the flame. Note that solid parts of the work to be soldered should be preheated first. During soldering, care must be taken to not remove the flame from the soldering area. The solder is to be introduced only when the solder point has been heated until it glows white-hot (Fig. 3). The actual soldering process should be carried out as quickly as possible. After the soldering process has been completed (Fig. 4), the bridge should cool down slowly.

It must in no case be quenched.

7. Contraindications and adverse reactions

Signs of intolerance to CoCrMo-Sold 1 are extremely rare if the manufacturer’s Instructions for use are adhered to. If the patient has a proven allergy against any component within the alloy, this alloy must not be used for safety reasons. There have been individual reportings of local irritations which were electrochemically induced. If various alloys have been used, it is possible that galvanic effects may occur.

Exposure to dust and vapors may cause irritation of the eyes and/or respiratory tract. Individual alloy components may have a carcinogenic effect. Please adhere to warnings and safety data sheets. In the State of California elements of this alloy are considered to be carcinogenic.

8. Warnings and precautions to be taken

Mechanical processing of the solder can lead to the formation of metal dusts. The mechanical process will also create heat and possibly burs. Processed material may therefore be hot and/or have sharp edges. Processing at temperatures above the solidus temperature may generate vapors. The workplace should be fitted with a suitable suction or ventilation system. We generally recommend that gloves, protective clothing, safety glasses and respiratory protection equipment are worn. Should there be contact with the skin, please wash the affected area with soap and water for a few minutes. Should the product come into contact with the eye, please rinse the eye with a gentle stream of water or saline solution for at least 15 minutes. Please consult a doctor if you are exposed to individual elements.

Please observe the relevant safety data sheet (see www.dentaurum.com).

Interproximal or occlusal contact with other metallic elements and the use of different types of alloys in the same oral cavity should be avoided to prevent galvanic reactions or local irritations in individual cases. The alloy should not be ground or polished in the patient’s mouth. There is no scientific evidence available on the safety or efficacy of treatment of pregnant women or nursing mothers or children.

9. Reuse

Dental solder that has already been processed must not be re-used.

10. Conditions for storing

Normal changes to the surrounding conditions (e.g. temperature, pressure, light) have no negative impact on the product characteristics of cobalt-based alloys.

11. Disposal

Please adhere strictly to the rules that apply for your region, bearing in mind the details outlined in the safety data sheets.

12. Quality information

Dentaurum ensures faultless quality of its products. These recommendations are based upon Dentaurum’s own experiences. The dental professional is solely responsible for the processing of the products. Responsibility for failures cannot be taken, as we, Dentaurum, have no influence on the processing on site.

13. Further information


Should the dental professional and/or the patient become aware of serious problems arising from the use of the product, it is important that the manufacturer and the competent authority in the country in which the dental professional and/or the patient is resident is informed accordingly.

The SSCP is available at https://ec.europa.eu/tools/eudamed.

14. Explanation of symbols used on the label

⚠ Please refer to the label. Additional information can be found at www.dentaurum.com (Explanation of symbols REF 989-313-00).

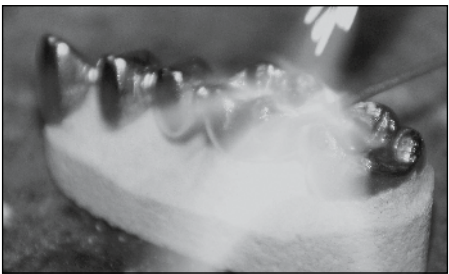Fig. 3

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