

Cleaning – so easy.


The Wash-Tray developed in collaboration with PROFESSIONAL.





Cleaning – so easy.

tioLogic® easyClean – the tray for machine processing.1

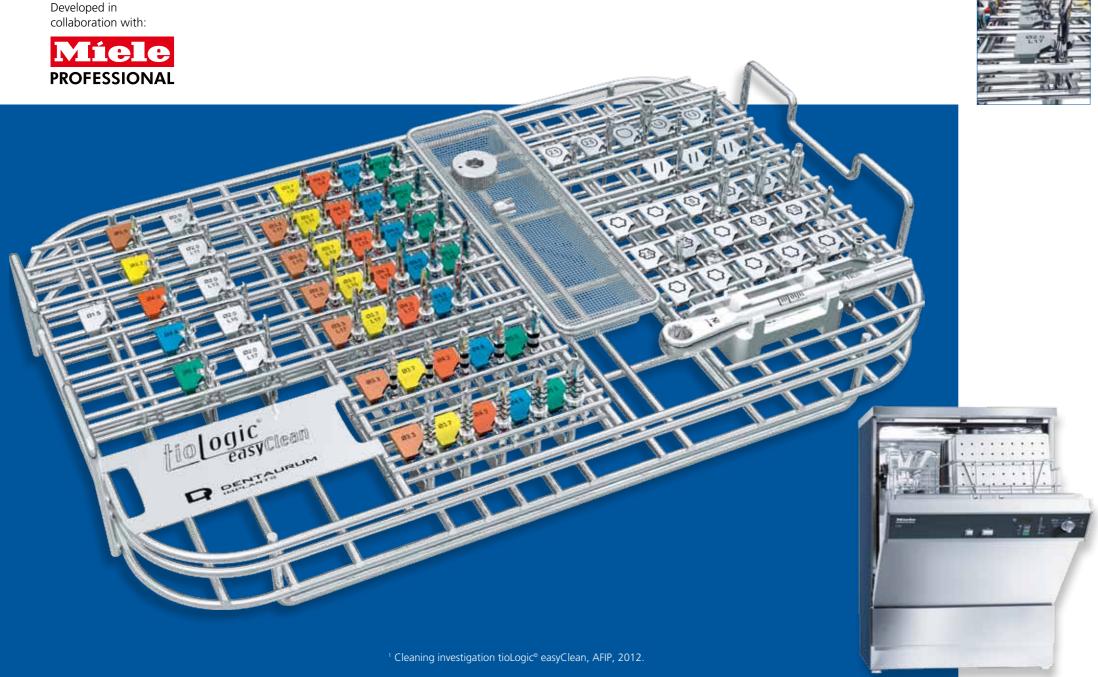

Innovative lattice structure



Stable fixation spring.



Reproducible machine-treatment guidance system. results.



The tioLogic® easyClean is a Wash-Tray that contains all the rotary instruments and accessory components needed for implantation. They are arranged according to the operation sequence. For optimal orientation there is a colour-coded and laser-labelled plastic clip beside each instrument. The used instruments and accessory components are put back in the corresponding slots directly after each use. This increases safety during implantation as all instruments are always located in the intended place. After implantation, the completely packed tioLogic® easyClean is transferred to the machine treatment cycle. Small parts and accessory components to be disassembled are placed in the mesh tray.

Machine treatment of the completely packed surgical tray.

Guaranteed by:

- optimal washing of the components
- minimal contact surface
- innovative lattice structure
- fine mesh tray

Added value for the practice:

- no manual treatment necessary
- reproducible machine-treatment results
- saves working time
- no sorting of components necessary

Suitable for doctors' surgeries and central treatment centres.

3



tioLogic® easyClean Wash-Tray

Fully equipped for placement of tioLogic® implants.

- for ø 3.3 mm / ø 3.7 mm / ø 4.2 mm / ø 4.8 mm / ø 5.5 mm
- reusable instruments
- standard aluminium container with tioLogic® easyClean Wash-Tray for cleaning and disinfecting machines
- hydrophob PTFE filter
- incl. tioLogic[®] DVD
- incl. tioLogic[®] Surgery Manual

tioLogic® easyClean (without contents)

REF 387-419-00

REF 387-418-00

Wash tray (without contents) for cleaning and disinfecting machines.



Further Information

Phone +497231/803-560 www.dentaurum-implants.de

